




INSTRUMENTS FOR CORNEAL PROCEDURES

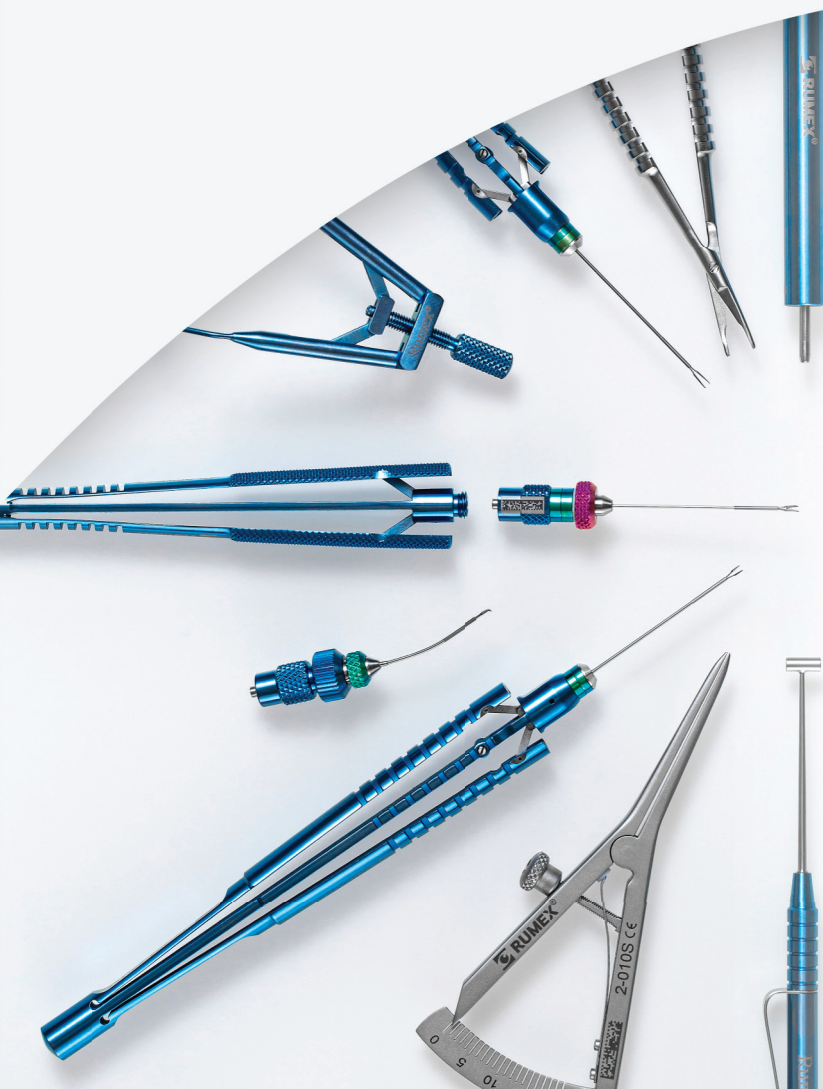

AFFORDABLE
INNOVATIONS

INTRODUCTION

The last decades have brought a revolutionary shift in the treatment of corneal endothelial disease. 20 years ago the only surgical treatment for corneal disorders was penetrating keratoplasty.

Although used successfully for over a century, penetrating keratoplasty requires many months of refractive adjustments before the eye achieves visual stability. Starting with the advent of posterior lamellar keratoplasty, a number of procedures have been developed, refined, and widely adopted, which have given patients faster recoveries and improved globe stability in comparison to traditional corneal transplantation. Each iteration of endothelial keratoplasty has involved the increasingly selective transplantation of corneal endothelial cells.

Constantly tracking the latest developments in eye surgery, we at RUMEX have designed a lineup of top quality instruments for all up to date techniques of corneal transplantation.

All instruments in this brochure are organized as sets according to the procedures of penetrating and lamellar keratoplasty. Sets are advisory, and if your surgical technique or professional preference requires a customized set, you can easily modify it by picking necessary instruments from a wide range of RUMEX products.

HOW TO PLACE AN ORDER

All the orders can be easily placed and paid online according to your shipping location directly via one of RUMEX online stores or through your local RUMEX representatives.



If your shipping country is USA you are welcome to order at **rumex.us**



If your shipping country is within Europe please visit **rumex.eu** to make your order.



For other regions visit **rumex.com** and find your local RUMEX representative.

Pricing

The prices are detected automatically by your location. The pricing policy may vary from region to region. If you are an international customer, please ask your local distributor for the current prices.

Shipping

We provide our retail customers with two delivery options: via local distributor or by direct shipment from our warehouses. Purchasing with our company is simple and convenient. Processing orders quickly and efficiently is a matter of primary importance to us!

Warranty conditions

For all instruments, RUMEX provides a lifetime warranty against any manufacturing or material defects. After carrying out a due expert analysis, if the defect was not caused by the improper handling or misuse, we will provide you either a 100% compensation or a free of charge exchange of a defective instrument for a new one. In some cases when instruments are improperly used or mishandled this may lead to occurrence of non-manufacturing defects which are not covered by RUMEX lifetime warranty. To avoid such cases please read carefully and always follow our sterilization and care instructions or consult our customer service for proper handling instructions.

+1 727 535 9600 (for USA, Canada)

+371 6616 3182 (for Europe, Asia, Africa, Latin America)

CONTENTS

› Cooperation with Eric Abdullayev, M.D., MBA, CEBT	4
› FEATURED	
› DALK	
› “Big Bubble” Creation	5
› Dissection and Separation	6
› Quadrisection	6
› Endothelial Transplantation (DSEK, DSAEK, DMEK)	
› Graft Preparation	7
› Marking	9
› Dissection	10
› Descemet’s Membrane Removal	10
› Implantation of the Graft	14
› Deep Lamellar Endothelial Keratoplasty (DLEK)	15
› BASIC	
› CORNEAL Transplantation	
› Markers	16
› Diamond Knives	17
› Corneal Trephine Blades	17
› Keratometer	17
› Scissors	18
› Forceps	19
› Miscellaneous	21
› SETS	
› DALK	22
› DLEK	23
› DSEK, DSAEK, DMEK	24
› Basic Corneal Instrument Set	25

INSTRUMENTS FOR CORNEAL PROCEDURES

COOPERATION WITH ERIC ABDULLAYEV, M.D., MBA, CEBT

Some of the featured corneal instruments were developed in cooperation with Eric Abdullaev, M.D., MBA, CEBT – Manager of Clinical Development and Innovations at Lions World Vision Institute, USA. He's been involved in eye banking community since 1996, after ten years of being cardiothoracic.

Dr. Abdullaev has many scientific publications on new techniques and procedures of eye banking, corneal transplantation and preservation. He is a creator of Nano-Cut Precise™ and Punch Ready DMEK™ ocular grafts.



3-0231



3-024T



4-261S



4-254S



3-0231	Abdullaev Corneal Marker for Keratoplasty Double-Ended (10.00 mm and 11.00 mm Diameters), with Central Marking Point, Round Titanium Handle
3-024T	Abdullaev I & II Marker (for DSAEK/DMEK Grafts) Titanium
4-261S	Abdullaev DMEK Grasping Forceps Stainless Steel
4-254S	Lambright-Abdullaev Ultrathin DSAEK Grasping/Inserting Forceps Stainless Steel


DEEP ANTERIOR LAMELLAR KERATOPLASTY (DALK)

Partial-thickness cornea transplant procedure allows to perform selective transplantation of the corneal stroma, leaving the native Descemet's membrane and endothelium in place.

“BIG BUBBLE” CREATION

13-172

Dissector for DALK Procedure

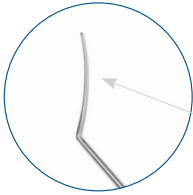


- Blunt beveled tip helps to create a track
- in deep stroma for the further cannula inserting.



15-450-27

Cannula for DALK Procedure, 27 Ga

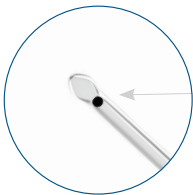

- Bottom port 0.20 mm
- Allows to achieve ideal “Big Bubble”




DEEP ANTERIOR LAMELLAR KERATOPLASTY (DALK)

DISSECTION AND SEPARATION

13-170



Trisector for DALK Procedure

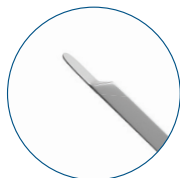

- Facilitates separation of stromal attachments from Descemet's membrane
- Flat tip completes dissection in the "Big Bubble" technique

13-171

Spatula for DALK Procedure

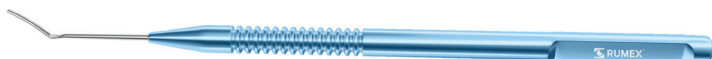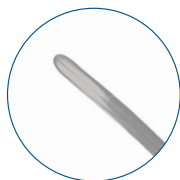

- Designed to finish incomplete dissection
- Center groove guides the blade for controlled stromal enlargement

QUADRISECTION

Matte finish/anti-glare coating

Minimal friction


Maraging stainless steel

Unique hardening system promoting 3,000+cuts without resharpening

Sharp cutting edges

11-038S Right
11-0381S Left (shown)


DALK Corneal Transplant Scissors

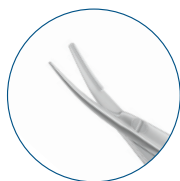

- Blunt ledge
- Protects Descemet's membrane

11-134 Right (shown)
11-135 Left

Holland Spatulated DALK Scissors

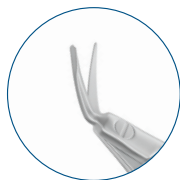

- Vannas-style micro blades with steep angle facilitate stromal dissection to the periphery
- Rounded blunt tips protect underlying tissue during cutting

DLEK, DSEK, DSAEK, DMEK

ENDOTHELIAL TRANSPLANTATION (DMEK, DSAEK, DSEK)

Endothelial transplantation is a surgical procedure that involves replacement of diseased posterior cornea with donor tissue while retaining the anterior corneal layers. The aim is to replace only the dysfunctional endothelial layer with healthy functioning endothelium.

DLEK (Deep Lamellar Endothelial Keratoplasty)

DSEK (Descemet's Stripping Endothelial Keratoplasty)

DSAEK (Descemet's Stripping Automated Endothelial Keratoplasty)

DMEK (Descemet's Membrane Endothelial Keratoplasty)

GRAFT PREPARATION


13-185

Tan Marginal DMEK Dissector



- Double-tipped end safely cuts peripheral donor Descemet's membrane (DM) without radial tears

- Curved single tip enables precise separation of DM from stroma

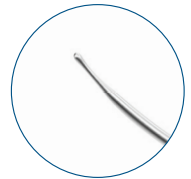

3-208T

DMEK/DSAEK "S" Marker



- Thin "S" leaves a fine line with minimal ink transfer to corneal stroma
- Optimal angle ensures better visualization during marking

DLEK, DSEK, DSAEK, DMEK

GRAFT PREPARATION

Designed in cooperation
with Eric Abdullayev, M.D., CEBT,
Lions Eye Institute for Transplant &
Research, Inc., USA.



3-0231 Abdullayev Corneal Marker for Keratoplasty




- Double-ended (10.00 mm and 11.00 mm diameters)
- Central marking point for precision
- Improves centration of the cornea during DSAEK microkeratome preparation
- Allows quick and accurate central dot placement
- Facilitates positioning of donor cornea on the donor punch

3-024T Abdullayev I & II Marker for DSAEK/DMEK Grafts



- I & II marks do not interfere with vision
- More stable staining
- Allows to apply 1.50 mm straight I & II marks at the very edge of the graft
- No additional tissue manipulation (no punch holes, no folding or unfolding of the graft)
- Saves time during graft preparation

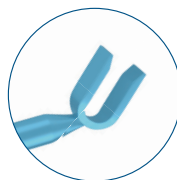


INSTRUMENTS FOR CORNEAL PROCEDURES

DLEK, DSEK, DSAEK, DMEK

4-240

Guell DMEK Forceps

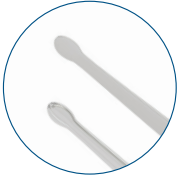

- Designed for safe removal of endothelium from donor cornea
- Highly polished broad tips allow peeling endothelium without tearing

4-261S

Abdullayev DMEK Grasping Forceps



Designed in cooperation
with Eric Abdullayev, M.D., CEBT,
Lions Eye Institute for Transplant &
Research, Inc., USA.



- Improved platform angle ensures relaxed hand position and better control
- Thin horizontal platforms provide stable membrane holding during separation

MARKING

3-204T

John DSAEK Double-Ended Marker 8.00/9.00mm



- Used for marking both recipient and donor cornea
- Combines two popular diameters: 8.00 mm and 9.00 mm
- Circular mark on recipient's cornea serves as guiding reference for Descemetorhexis


INSTRUMENTS FOR CORNEAL PROCEDURES

DLEK, DSEK, DSAEK, DMEK

DISSECTION

13-137 Straight (shown)
13-138 Curved

Corneal Dissector

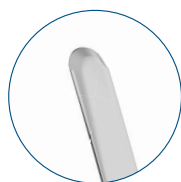

Straight

- For efficient intrastromal dissection



Curved


- Designed for proximal corneal stroma dissection
- Blade contour follows corneal curvature for precise stromal work



DESCEMET'S MEMBRANE REMOVAL

5-0322

Reversed Sinskey Hook

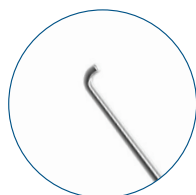

0.15 mm tip

- Ideal for scoring the recipient bed and placing donor lamella
- Blunt tip breaks through Descemet's membrane gently and safely
- Provides soft membrane segmentation without damage
- Ensures easy localization of donor's lamella
- 0.15 mm tip for fine access



13-139/I

Irrigating Endothelial Stripper



- Thin, semi-circular tip angled 90° from the irrigation tube
- Used to gently peel and remove the endothelial layer




DLEK, DSEK, DSAEK, DMEK

DESCEMET'S MEMBRANE REMOVAL


The DMEK forceps have been designed to assist in DMEK procedures. The forceps are used to peel off the Descemet's membrane once it has been scored with the Descemet's spatula.

4-251

Cross-Action DMEK Forceps

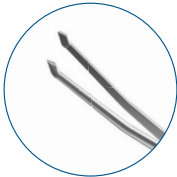


- Cross-action prevents hyperextension of incision
- Reversed delicate tips with tungsten carbide for enhanced grip

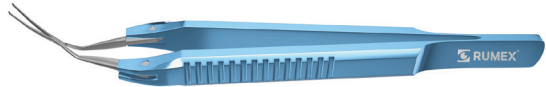


4-262

AVAILABLE
MODIFICATION

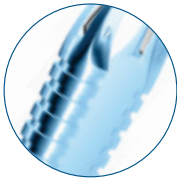


- Lightweight and ergonomic round titanium handle option

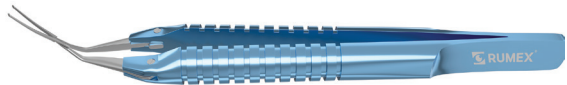

4-246S

Florakis Endothelial Forceps

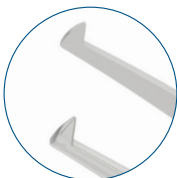

- Reversed triangular tips angled at 75° for precision and visualization
- Tungsten carbide coated platform ensures secure grip



4-247*

Florakis Microinvasive Endothelial Forceps



- Designed to separate endothelium from recipient cornea
- Reversed triangular tips for controlled handling
- 23 Ga tube enables procedures through ≤ 1.00 mm incision
- Allows access to all corneal areas without hyperextending the incision



*Tip Only. To be used with Universal Handle 12-003T

Product design and/or features that do not influence its functionality and main parameters are subject to change.

DLEK, DSEK, DSAEK, DMEK

DESCEMET'S MEMBRANE REMOVAL

Specific design of the spatula contributes to a free access to almost all parts of the recipient's inner cornea.



13-1491


John DSAEK Descemet's Stripper



- Provides excellent contact with the inner cornea
- Ideal for cases with strong adhesion of Descemet's membrane to the recipient's corneal stroma

13-182

John Dexatome DMEK/DSAEK Spatula

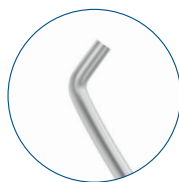

- Designed for easy removal of Descemet's membrane as a single disc

13-183

John DSAEK Stromal Scrubber



- Sandblasted tip roughens the stroma in peripheral areas of the Descemetorhexis
- Enhances graft adherence to the patient's cornea
- Minimizes the risk of graft detachment

13-184

John DSAEK Glider



- Facilitates smooth gliding of the donor disc
- Clears fluid and smoothes the corneal surface at the donor-recipient interface

DLEK, DSEK, DSAEK, DMEK

DESCEMET'S MEMBRANE REMOVAL

13-151S

Cindy Sweeper DSEK Spatula



- Smoothens the area between donor graft and recipient's stromal bed
- Gently glides over the recipient's cornea

13-152S

Carlson DSEK Smoother



- Removes residual fluid between donor graft and recipient's stromal bed
- Helps secure graft position post-implantation by gentle surface pressure

13-153S

Terry DSEK Scraper




- Roughens the inner corneal stroma of the recipient
- Promotes better graft adherence

**13-154T
13-155T**

Melles DSAEK PLK Scraper




For efficient stripping of the Descemet's membrane

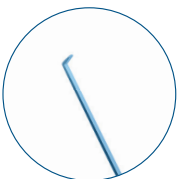

45° angled shaft
45° angled tip

13-154T

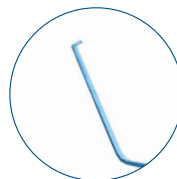

45° angled shaft
90° angled tip


13-155T

DLEK, DSEK, DSAEK, DMEK

IMPLANTATION OF THE GRAFT

13-150T


Spatula-Guide for Corneal Endothelium Implantation

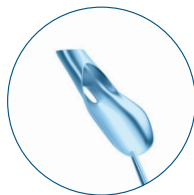

- Prepares the graft for insertion while protecting endothelial cells during implantation
- Prevents compression of the donor disc as it passes into the inner cornea
- Can be used together with 4-034 Forceps for Corneal Endothelium Implantation

4-034*



Forceps for Corneal Endothelium Implantation

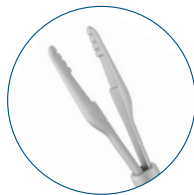

- Designed for inserting the donor button using the pull-through technique
- Ensures secure placement in the stromal layer, minimizing endothelial cell damage
- Blunt tips with delicate atraumatic grooves for firm fixation of the corneal stromal layer

4-2019T

Corneal Donor Insertion Forceps

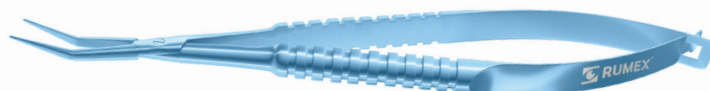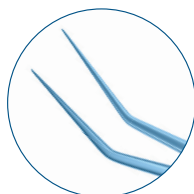

- Designed for atraumatic insertion of the donor lamella folded in a taco shape
- 30° angled, 18.00 mm tip for optimal access and control

*Tip Only. To be used with Universal Handle 12-003T

DLEK, DSEK, DSAEK, DMEK

IMPLANTATION OF THE GRAFT

Designed in cooperation with
Eric Abdullayev, M.D., CEBT,
Lions Eye Institute for Transplant &
Research, Inc., USA.




4-254S


Lambright-Abdullayev Ultrathin DSAEK Grasping/Inserting Forceps



Optimized for inserting ultrathin DSAEK grafts (thickness 70 microns or less)

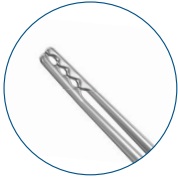

- Optimized for inserting ultrathin DSAEK grafts (thickness 70 microns or less)

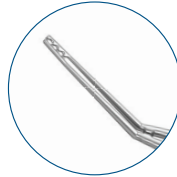

- Wave-shaped serrations ensure non-slip insertion and reduce risk of tissue adherence

DEEP LAMELLAR ENDOTHELIAL KERATOPLASTY (DLEK)

Matte finish/anti-glare coating

Minimal friction

Maraging stainless steel

Unique hardening system promoting 3,000+cuts without resharpening

Sharp cutting edges

DLEK Scissors

Designed for removal
of the posterior cornea



11-036S

Medium curve



11-0361S


Strong curve

CORNEAL TRANSPLANTATION


MARKERS

3-140T

Corneal Transplant Marker

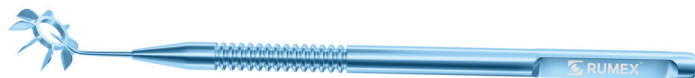


- 7.00 mm ring with 8 radial blades
- Center pointer for improved centration without obstructing visibility

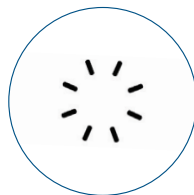


3-0304T

Osher-Neumann Corneal Marker

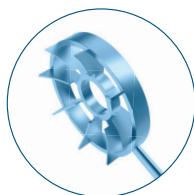

- Low-profile design
- Equipped with 8 radial blades

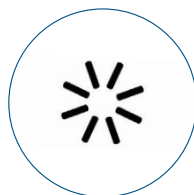


3-0219T

Hoffer Optical Zone Marker



- Low-profile design
- Equipped with 8 radial blades

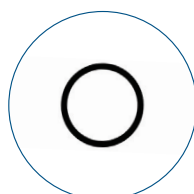

CORNEAL TRANSPLANTATION

DIAMOND KNIVES

Extreme sharpness


Minimal pressure

Reduced risk of wound leak

Predictable cuts and architecture of the wound for more consistent healing

6-10/6-053


Side-Port Diamond Knife, Trifacet Blade

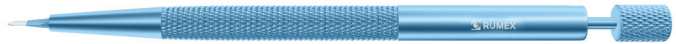

- Trifacet blade for stab incisions, initial groove creation, and tunneling
- Blade width: 1.00 mm

6-20/6-107

Phaco Diamond Knife, Trapezoid Self-Diving Blade

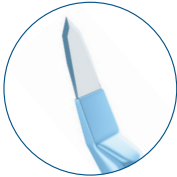

- Symmetrical sharp edges for precise, distortion-free incisions
- Enables smooth self-diving into tissue
- Blade width: 2.00 / 2.30 mm

CORNEAL TREPHINE BLADES

16-0300 6.00 mm

16-0307 8.00 mm

16-0301 6.50 mm

16-0308 8.25 mm

16-0303 7.00 mm

16-0309 8.50 mm

16-0305 7.50 mm

16-0310 9.00 mm

16-0306 7.75 mm


16-0311 9.50 mm



16-020T

Maloney Intraoperative Keratometer

Designed for qualitative assessment of astigmatism after donor cornea suturing

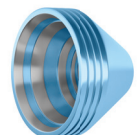

CORNEAL TRANSPLANTATION

SCISSORS

Matte finish/anti-glare coating

Minimal friction

Maraging stainless steel

Unique hardening system promoting 3,000+cuts without resharpening

Sharp cutting edges

11-040S

Westcott Tenotomy Scissors



- Blunt atraumatic tips



11-044S

Westcott Stitch Scissors



- Pointed tips



11-020S Left
11-0201S Right (shown)

Katzin Corneal Transplant Scissors




- Strongly curved delicate blades
- Designed for precise cutting of corneal tissue



11-024S Left (shown)
11-0241S Right

Castroviejo Corneal Section Scissors

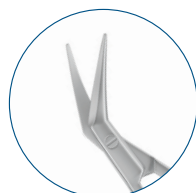

- Lower blade is 0.50 mm longer than the upper blade
- Allows minimal corneal alteration, adapting to its curvature and thickness




CORNEAL TRANSPLANTATION

FORCEPS

5.00 mm tying platform hard-faced with tungsten carbide coating for even greater grip and increased control

4-0505T Micro Colibri Corneal Forceps

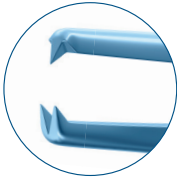

- 1×2 teeth, 0.12 mm




4-0541T Castroviejo Colibri Corneal Forceps




- 1×2 teeth, 0.12 mm



4-0814T Pollack Corneal Transplantation Fixation Forceps

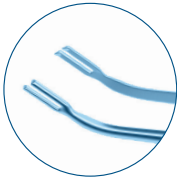

- Pierce-type tips for two-point fixation
- Double prongs securely hold the cornea
- Delicate 1×2 teeth, 0.12 mm
- 1.2 mm spacing between tips

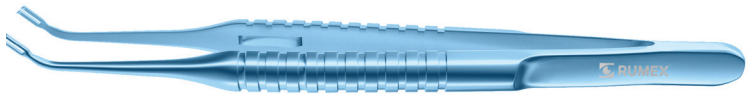


CORNEAL TRANSPLANTATION

FORCEPS

Tying platform hard-faced with tungsten carbide coating for even greater grip and increased control over the tissue and suturing materials

4-0600S
4-0600T


Castroviejo Suturing Forceps

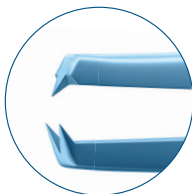

- 1×2 teeth, 0.12 mm tip
- Tying platform, 6.00 mm

4-0607S

Bishop-Harmon Suturing Forceps


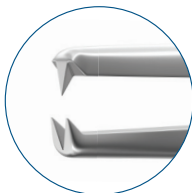


- 1×2 teeth, 0.30 mm tip
- Tying platform, 5.00 mm

4-178S
4-174T

McPherson Tying Forceps

FOR 8.0 TO 11.0 SUTURES

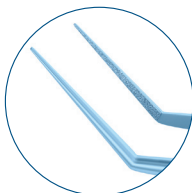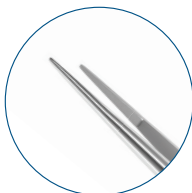

- Angled shafts allow smooth maneuvering in the anterior chamber
 - Tying platform, 8.00 mm
 - Titanium
- 4-174T**

CORNEAL TRANSPLANTATION

MISCELLANEOUS


Tungsten carbide coated tips for advanced gripping

Lightweight

Non-magnetic

8-031T
8-045T


Barraquer Needle Holder

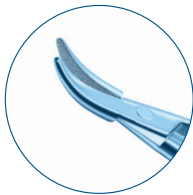

FOR 4.0 TO 7.0 SUTURES

- Standard jaws 12.00 mm, without lock

8-031T

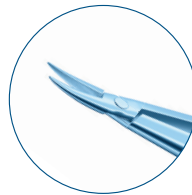

FOR 8.0 TO 11.0 SUTURES

- Extra fine jaws 8.00 mm, without lock

8-045T

13-110


Paton Double-Ended Spatula and Spoon



- Teardrop-shaped spoon ensures safe corneal button transfer
- Spatula promotes delicate epithelium manipulations

15-051-25

Rycroft Anterior Chamber Cannula, 25 Ga

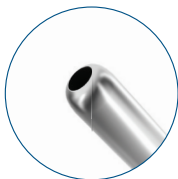

Smooth blunt tip



- Used to maintain and form the anterior chamber by injecting or removing air, fluids, viscoelastics, and intraocular medications
- Smooth blunt tip
- 4.00 mm length, 45° angle

15-301/303


Silicone Bulb with Adapter



INSTRUMENTS FOR CORNEAL PROCEDURES

FEATURED SETS

DALK SET

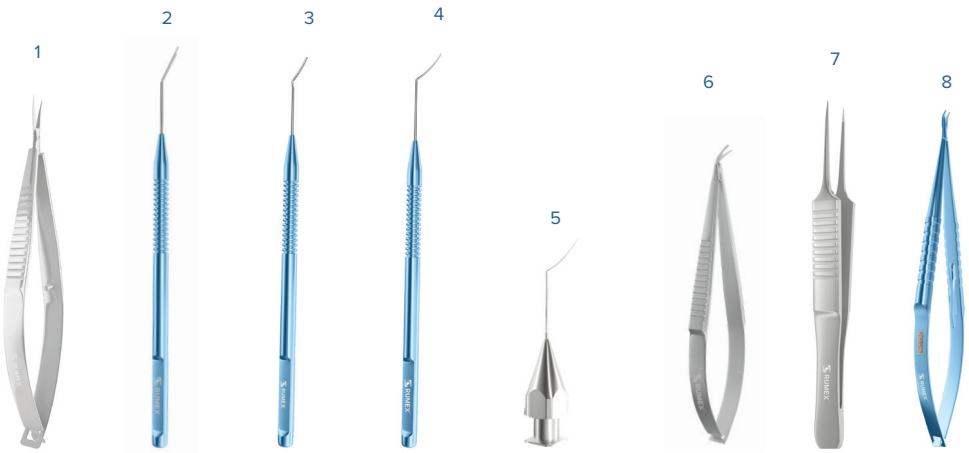

Key	Reference	Description
1	11-038S 11-0381S*	Scissors for DALK Procedure, Right Scissors for DALK Procedure, Left
2	13-170	Trisector for DALK Procedure
3	13-171	Spatula for DALK Procedure
4	13-172	Dissector for DALK Procedure
5	15-450-27	Cannula for DALK Procedure, 27 Ga
6	11-134	Holland Spatulated DALK Scissors, Right
	11-135*	Holland Spatulated DALK Scissors, Left
7	4-178S	McPherson Tying Forceps, Straight
	4-174T*	McPherson Tying Forceps, Angled
8	8-024T	Barraquer Needle Holder, Extra Fine Jaws, with Lock


*not shown



Plastic Sterilization Tray


with finger tip mat
254×152×38 mm
10×6×1.5 in

18-305

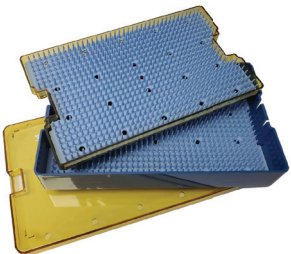

FEATURED SETS

DLEK SET

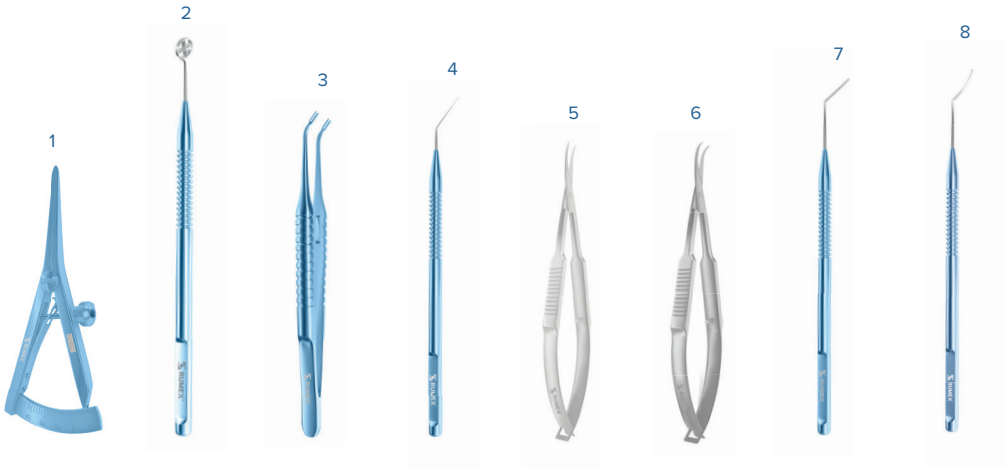


Key	Reference	Description
1	2-010T	Castroviejo Caliper
2	3-0217T	Hoffer Optical Zone Marker, 8.00 mm
3	4-0814T	Pollack Corneal Transplantation Fixation Forceps
4	5-0322	Reversed Sinskey Hook
5	11-036S	DLEK Scissors, Medium Curve
6	11-0361S	DLEK Scissors, Strong Curve
7	13-137	Corneal Dissector, Straight
8	13-138	Corneal Dissector, Curved



Plastic Sterilization Tray


with finger tip mat
254x152x19 mm
10×6×0.75 in

18-304

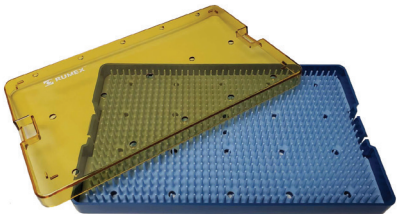

FEATURED SETS

DSEK, DSAEK, DMEK SET

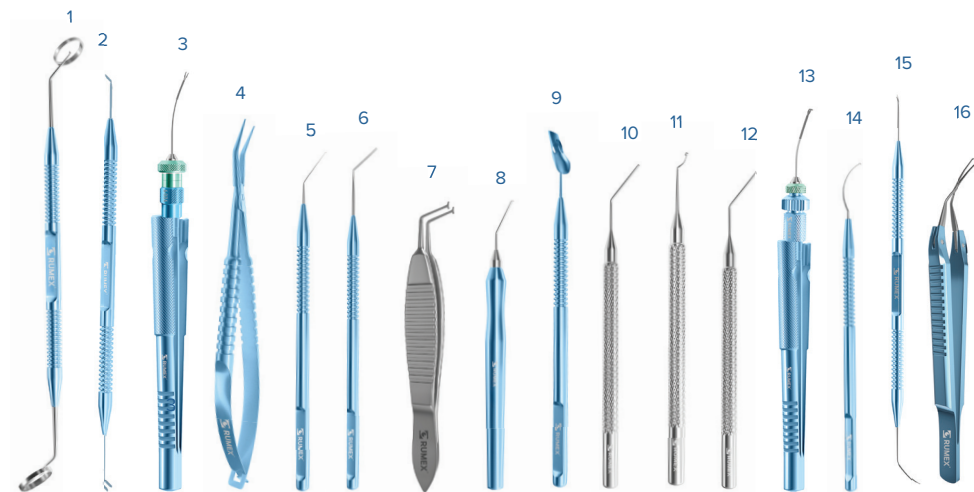


Key	Reference	Description
1	3-0231	Abdullayev Corneal Marker for Keratoplasty
2	3-024T	Abdullayev I & II Marker (for DSAEK/DMEK Grafts)
	3-208T*	DMEK/DSAEK "S" Marker
3	4-034	Forceps for Corneal Endothelium Implantation, 23 Ga/
	/12-003T	Universal Squeeze Handle
4	4-2019T	Corneal Donor Insertion Forceps
5	5-0322	Reversed Sinskey Hook
6	13-137	Corneal Dissector, Straight
	13-138*	Corneal Dissector, Curved
7	4-261S	Abdullayev DMEK Grasping Forceps
8	13-139/I	Endothelial Stripper, Irrigating
9	13-150T	Spatula-Guide for Corneal Endothelium Implantation
10	13-151S	Cindy Sweeper DSEK Spatula
11	13-152S	Carlson DSEK Smoother
12	13-153S	Terry DSEK Scraper
13	4-247	Florakis Microinvasive Endothelial Forceps, 23 Ga/
	/12-003T	Universal Squeeze Handle
	4-246S*	Florakis Endothelial Forceps
14	13-182	John Dextatome DMEK/DSAEK Spatula
	13-183*	John DSAEK Stromal Scrubber
	13-184*	John DSAEK Glider
15	13-185	Tan Marginal DMEK Dissector
16	4-251	Cross-Action DMEK Forceps
	18-305*	Plastic Sterilization Tray with Silicone Finger Mat, Double Level, Extra Large

*not shown

FEATURED SETS

CORNEAL TRANSPLANTATION SET

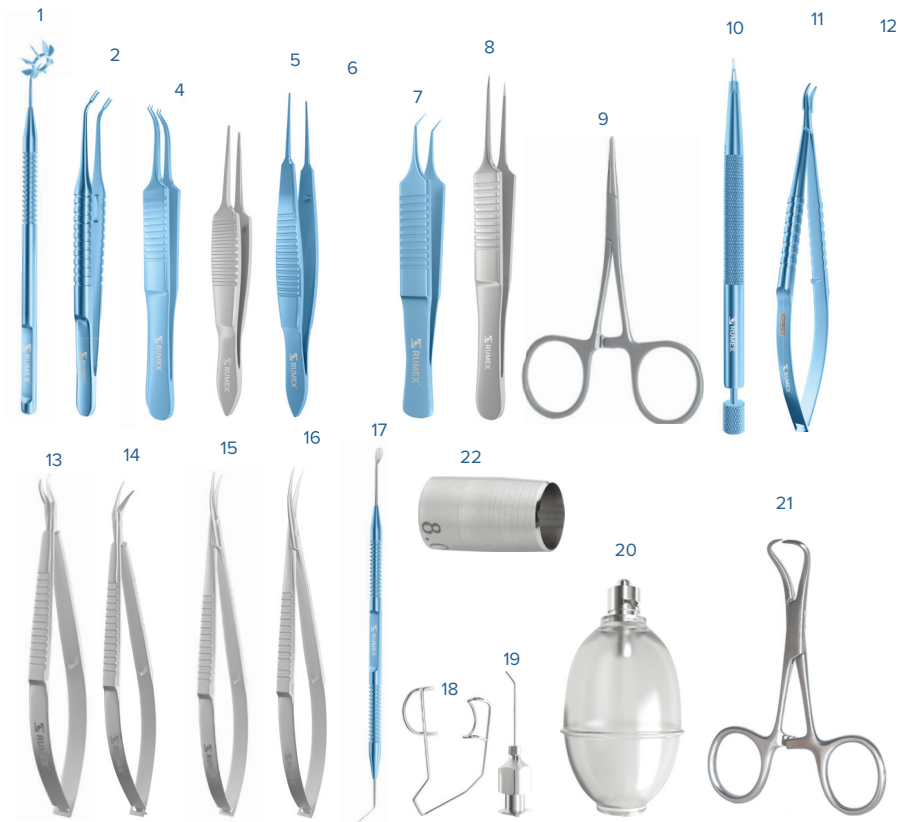

Key	Reference	Description	Key	Reference	Description
1	3-140T	Corneal Transplant Marker	11	8-031T	Barraquer Needle Holder, Standard Jaws, without Lock
2	4-0814T	Pollack Corneal Transplantation Fixation Forceps	12	8-045T	Barraquer Needle Holder, Extra Fine Jaws, without Lock
3	4-0541T	Castroviejo Colibri Corneal Forceps, 0.12 mm, 1x2 Teeth	13	11-020S	Katzin Corneal Transplant Scissors, Left
4	4-0607S	Bishop-Harmon Suturing Forceps, 0.30 mm, 1x2 Teeth	14	11-0201S*	Katzin Corneal Transplant Scissors, Right
5	4-0600T	Castroviejo Suturing Forceps, 0.12 mm, 1x2 Teeth	15	11-024S	Castroviejo Corneal Section Scissors, Right
6	4-0601T	Castroviejo Suturing Forceps, 0.30 mm, 1x2 Teeth	16	11-0241S*	Castroviejo Corneal Section Scissors, Left
7	4-090T	Kelman-McPherson Tying Forceps, 4 mm Platforms	17	11-040S	Westcott Curved Tenotomy Scissors
8	4-178S	McPherson Straight Tying Forceps, 7 mm Platforms	18	11-044S	Westcott Curved Tenotomy Scissors
9	4-120S	Hartman Mosquito Forceps	19	13-110	Paton Double Ended Spatula And Spoon
10	6-10/6-053	Diamond Knife, Trifacet Blade, 1.00 mm	20	14-022S	Barraquer Wire Speculum, Adult Size
	6-20/6-107*	Diamond Knife, Self-Diving Trapezoid Blade, 2.00/2.30	21	15-051-25	Rycroft Anterior Chamber Cannula, 25 Ga
			22	15-301/303	Silicone Bulb With Adapter
				16-180S	Towel Forceps
				16-0305*	Corneal Trephine Blades, 7.50 mm
				16-0307	Corneal Trephine Blades, 8.00 mm
				18-305*	Plastic Sterilizing Tray, Double Level, Large

*not shown



RUMEX INTERNATIONAL COMPANY

 **14240 Carlson Circle,
Building K, Suite 8, Tampa,
FL 33626, USA**

RUMEX BALTICS SIA

 **Robezu iela 46, Riga,
Latvia, LV-1004**

USA & CANADA:

 **+1 (727) 535 9600**
 **+1 (877) 77 RUMEX (Toll-Free)**
 **+1 (727) 535 8300**
usacs@rumex.com

EUROPE, ASIA, AFRICA, LATIN AMERICA:

 **+371 6616 3182**
 **customerservice@rumex.com**

 **rumex@rumex.com**
 **www.rumex.com**

